




CUSTOMER: KURITA EUROPE GMBH, GERMANY

Kurita Europe is a subsidiary of Kurita Water Industries Ltd. headquartered in Tokyo, Japan. Founded in 1949, the Group has 6,600 employees worldwide with a turnover of approximately 2.3 billion USD. Kurita is one of the international market leaders in industrial water and process treatment. Founded 1989, the European team consists now of about 550 employees at 12 various locations, covering the EMEA (Europe, Middle East, Africa) sales area. With own production facilities Kurita contributes to a reliable supply chain for our customers. Own R&D centers, laboratories and pilot test plants allow Kurita to develop continuously customised state-of-the-art technologies.

www.kurita.eu

„The complete PTFE lining of the reactor as well as of the sensors for temperature and pressure measurement in combination with the simultaneous stirring have been a crucial factor for our purchase decision of the DB-300 reactor.“

- André Schirmer (R&D Specialist, Kurita Europe GmbH)

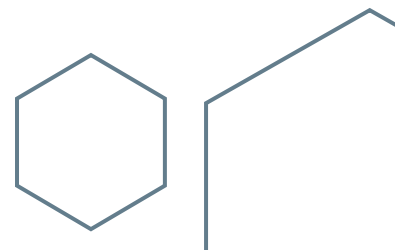



PRODUCTS + INSTRUMENTS

CHALLENGE

Due to the resistance of standard metal alloy, the researches on chemical and hydrodynamic processes in solutions with a high salt content in the field of energy technology and environmental technology have been a difficult problem to solve for Kurita Europe GmbH. The reactors also have to withstand increased temperatures and pressures for hours or even for days. In order to carry out comparative test series, alternating pressure load has been an issue effecting the resistance of the test apparatus, too.

APPLICATION

Kurita Europe GmbH uses the DB-300 reactor for tests in a temperature range between 100 and 220 °C. In addition, synthetically mixed salt solutions or salt solutions from customer specific processes are used. The researches mainly relate to product stability in the medium used as well as to corrosion tests based on metal coupons. Therefore, maintaining constant pressure and constant temperature conditions over several days is mandatory. For this process, the reactor should be easy to operate.

EQUIPMENT

highreactor **DB-300** reactor with BLH-650/BAH-300 combination for heating and stirring

RESULTS/BENEFITS

The complete PTFE lining of the reactor and the measuring sensors enable Kurita to work with strongly saline media. The continuous monitoring of pressure and temperature contributes to achieve high reproducibility in comparative tests. The design and the closure system of the autoclave allow an easy handling and impresses by its durability.

„The combination of the single components for different experimental setups allows us a multiple deployment of the DB-300 reactor and therefore a cost reduction. The easy handling of the often heavy reactors was another reason for the purchase decision”, explains André Schirmer, R&D Specialist at the Kurita Europe GmbH.

